



## Micro-Tec SH6 holder for single edge cutting blades

#52-004162

### Product description

The Micro-Tec SH6 holder for single edge cutting blades has been designed to hold single edge cutting blades with a width of 39mm. It has a large wooden handle and allows to use the blades for cutting or scraping applications. The holder clamps the blade and allows for safer handling of the single edge cutting blades.

The Micro-Tec SH6 blade holder is compatible with the following blades:

- #52-004150 Micro-Tec CB-S standard sharp carbon steel single edge cutting blade
- #52-004153 Micro-Tec CB-S standard sharp carbon steel single edge cutting blade in dispenser box
- #52-004155 Micro-Tec CB-D durable single edge carbon steel cutting blade
- #52-004160 Micro-Tec CB-H heavy duty single edge carbon steel cutting blade

All other single edge cutting blades with a width of 39mm and similar design as the cutting blades listed.

The Micro-Tec SH-6 blade holder is **not** suitable for holding razor blades.

### Directions for use

Use caution when changing blades – the blades do have a very sharp cutting edge.

Only hold the blades at the back or at the sides – do not touch the sharp edge.


Swing blade holder fully open to release the blade or to load a blade.

Locate the back of the blade in the groove

Close the blade holder completely to secure blade for use.

Reverse blade for safe storage.

No additional tools required for changing blades.

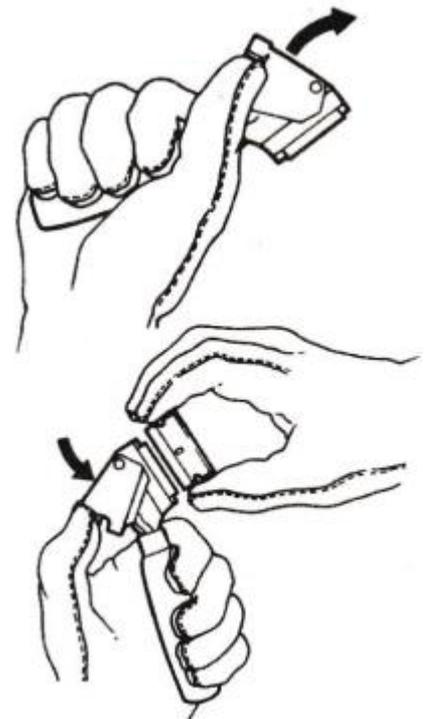

TSB 52-004162 Micro-Tec SH6 holder for single edge cutting blades, 2022-05-03 Revision 2

